

Sileo

Aged Care Bed

for Long Term
Care



A fully featured, easy to operate, electric bed designed for residential care facilities. With high-low range, safe and sturdy scissor lift mechanism, inbuilt base extension, removable head and foot panels and multiple positioning, this bed is designed for optimum support, strength and safety.

Available in Standard and King Single platforms.



novis
SPECIALIST BEDDING
SOLUTIONS

Sileo


Aged Care Bed

Sturdy, Slatted Frame

One piece mainframe constructed from powder coated steel for maximum strength. Operates as a high-low, four part profiling bed, with electric adjustable bed height, back rest, knee break and heel lift. Powered by a reliable Dewert actuator, with an inbuilt Back Up Battery available as an optional accessory.

The slatted design of the platform promotes the circulation of air in and around the mattress, helping to prevent moisture build up and regulate patient temperature – critical for both comfort and pressure injury prevention.

Available in either Standard or King single platforms.

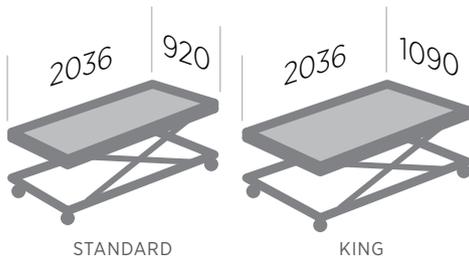

Scissor Lift Mechanism

The high-low scissor lift mechanism, with cross support beneath the platform, ensures the platform moves straight up and down, reducing the risk of damage to walls or other equipment around the bed. The design helps to evenly distribute the load, increasing durability and safety.



Central Locking Castors

Fitted with a central locking foot bar at the base of the bed, these smooth glide semi-enclosed castors are designed to minimise the build up of dirt and dust for easy cleaning and effortless operation.

Quick Release Panels

Both head and foot panels are easily removable, ideal for caregiver positioning when administering foot or hair care routines.

Three design options available. Classic curved design in two timber finishes or moulded polymer design with inbuilt handles and wall protectors.



Base Extension

The built in base extension increases the length of the bed platform by a maximum of 200 mm. The base easily extends at the foot end by loosening the extension knob and pulling the extension frame out to the desired length.




Hand Controller

Lightweight, water resistant hand controller with soft touch button for raising, lowering and positioning, including one touch auto contouring. Convenient inbuilt hanging hook for easy storage and access.



Electronic lockout function allows hand controller to be locked, ideal for patients incapable of operating the bed due to cognitive impairment. The lock is operated by a tamper proof magnetic swipe key attached to the base of the handset – simply swipe the key across the lock / unlock symbol to activate or deactivate the locking mechanism.

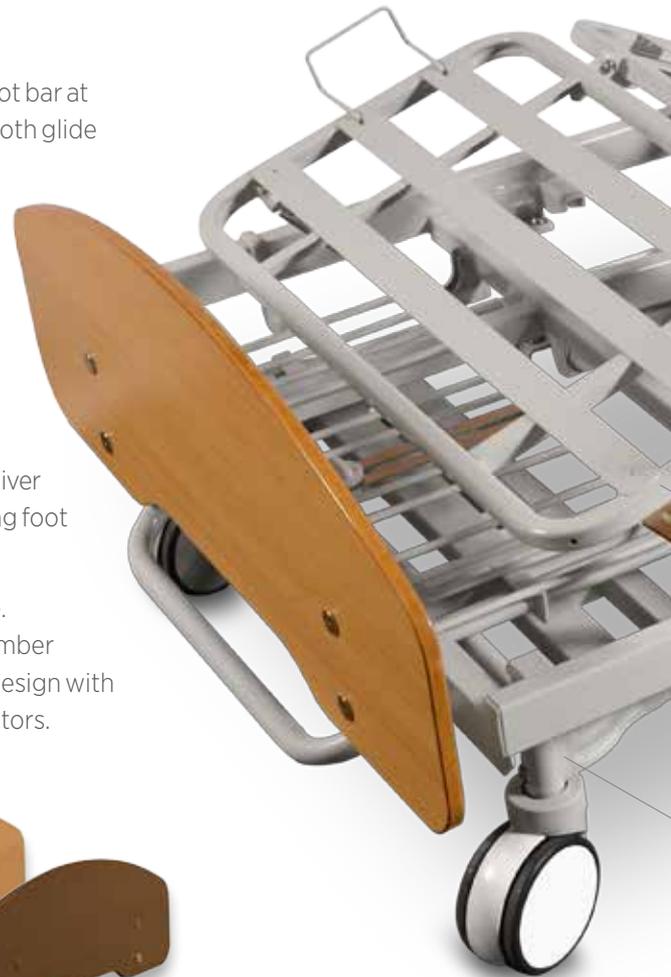



Four part profiling with Auto Contour

Drainage bag holders

Electronic hand controller with lockout

Detachable mattress retainers

Central locking foot bar

Accessories

Range of optional accessories available



Collapsible Side Rails

Collapsible metal side rail system attaches firmly to bed frame. Simply release the locking knob to fold down each rail for patient access or transfer.

CODE **BE6001** SOLD AS A PAIR



Back Up Battery

A rechargeable battery pack designed for a single operation, in the event of a power failure. Battery attaches to a cradle underneath the bed.

CODE **BE5001**



Self Help Pole

Support pole to aid patient repositioning while lying in bed. Mounts into inbuilt brackets at either side of bed head.

CODE **BE7001**

Sileo

Aged Care Bed



	SAFE WORKING LOAD	220 kg	
	MAXIMUM PATIENT WEIGHT	185 kg	
	WARRANTY	3 year Electrical / 10 year Frame / 1 year Back Up Battery	
	FRAME MATERIAL	Powder Coated Steel	
SYSTEM	MODEL	Standard	King
	OVERALL DIMENSIONS	W 920 mm [†] L 2036 mm [‡]	W 1090 mm [†] L 2036 mm [‡]
	HEIGHT BASE	330 - 780 mm	330 - 780 mm
	MATTRESS DIMENSIONS <small>RECOMMENDED</small>	W 860 - 900 mm L 1960 - 2200* mm	W 1050 - 1070 mm L 1960 - 2200* mm
	BED FRAME WEIGHT <small>WITHOUT ACCESSORIES</small>	102 kg	107 kg
CODE	BED WITH GOLDEN OAK PANELS	BE1001	BE2001
	BED WITH WARM OAK PANELS	BE1002	BE2002
	BED WITH MOULDED PANELS	BE1003	NA
COMPLIANCES	THERAPEUTIC GOODS ADMINISTRATION	ARTG 278782	ARTG 278802
	INTERNATIONAL STANDARDS ORGANISATION	ISO 9001	ISO 14001
	INTERNATIONAL ELECTROTECHNICAL COMMISSION	IEC 60601-2-52:2009	



[†] ± 5 mm
[‡] ± 10 mm
 * with bed fully extended

THERAPEUTIC DEVICES SHOULD ONLY BE USED IN ACCORDANCE WITH MANUFACTURERS INSTRUCTIONS AND UNDER THE CONSENT, SUPERVISION AND MANAGEMENT OF A SUITABLY QUALIFIED HEALTH PROFESSIONAL.

INFORMATION CORRECT AT TIME OF PUBLICATION - OCTOBER 2016 - © NOVIS HEALTHCARE, ALL RIGHTS RESERVED